





Environmental Triaxial Automated System (ETAS)

Overview: The Environmental Triaxial Automated System (ETAS) is a temperature controlled load frame-based triaxial testing system. Key features include:-

- Frozen Soil Testing
- · Gas Hydrate Testing
- · High Pressure Testing
- · High & Low Temperature Testing

The cooling system provides ideal conditions for frozen soil testing. Its high pressure testing capability (up to 100MPa) along with low temperature provides the ideal environment for Gas Hydrate Testing. The system also has the option for a Heating only system.

Key Features: Benefits to the User:

Heating & Cooling Options: -20°C to +65°C -30°C to +85°C -30°C to +100°C	The cooling system uses a coiled tube section inside the cell, this is connected through the standard ports within the cell base to a glycol cooling unit allowing temperature transfer close to the specimen.
Optional Heating Only System: Ambient to 60°C Ambient to 100°C	The heating system works using thermal pads attached to the outside of the triaxial cell then enclosed within an environmental chamber to retain the heat. There are up to 4 temperature sensors relaying back to the control box for enhanced temperature accuracy.
High Pressure (Up to 100MPa):	High Pressure combined with low temperature make it ideal for Gas Hydrate testing.
Gas Hydrate Testing (Gaseous back pressures):	The ETAS can be adapted for gaseous back pressures such as Air, Carbon Dioxide, Nitrogen or Methane for Gas Hydrate testing.
Frozen Soil Testing:	Cooling to -30°C provide ideal conditions for frozen soil testing.

Tests that can be Performed:

B-check, consolidated drained (CD) triaxial, consolidated undrained (CU) triaxial, consolidation (Triaxial), constant rate of loading (CRL), constant rate of strain (CRS), slow cyclic testing, K0 (K-Zero), multi-stage testing, quasi-static (low speed/creep) tests, stress paths and unconsolidated undrained (UU) triaxial. All tests can be performed with complete control of specimen temperature.

Upgrade Options:

Bender elements, unsaturated testing, local strain transducers and gaseous back pressures.

Technical Specification:

Load Range (kN):	50, 100, 250, 400, 500, 1000, 2000		
Pressure Range (MPa):	8, 14, 16, 20, 32, 64, 100		
Sample Sizes (mm):	38 to 150		
Available Temperature Ranges:			
Heating Only:	Ambient 65°C, Ambient to 100°C		
Heating & Cooling:	-20°C to +65°C, -30°C to +100°C, -30°C to +100°C		



Fig 1. System set-up with Gas Hydrate option

GDSLAB Software

GDSLAB is the control and data acquisition software. The software starts with a core application known as the kernel which, allows for data acquisition from your hardware, but no test control. Simply add the appropriate module or modules to complete the test suite functionality you require.

8 Channel Serial Pad


The standard GDS 8 channel data acquisition device, known as the "serial data pad", may be used within any system combination. This 16 bit device has 8 computer controlled gain ranges, designed to suit transducers used in a triaxial test. i.e.

- +/-10mV, +/-20mV, +/-30mV (load cells)
- +/-100mV, +/-200mV (pressure transducers)
- +/- 1V, +/- 5V, +/- 10V (displacement transducers)

Gaseous Back Pressures

- Carbon Dioxide
- Nitrogen
- Methane
- Air

Used in conjunction with a GDS PC Controlled pressure for gaseous pressure of up to 25MPa.

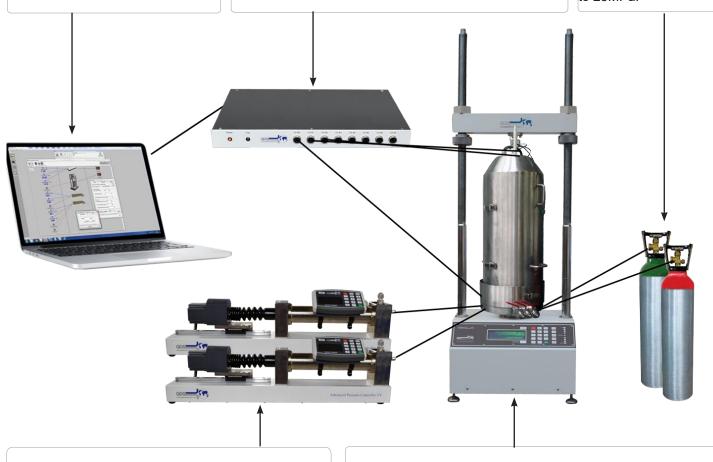

Pressure Volume / Controllers


- Enterprise Triaxial Automated System (ELETAS) which is based on 1MPa Enterprise Pressure / Volume Controller (ELDPC).
- Standard Triaxial Automated System (STDETAS) which is based on 3MPa Standard Pressure / Volume Controller (STDDPC).
- Advanced Triaxial Automated System (ADVETAS) which is based on 2MPa Advanced Pressure /Volume Controller (ADVDPC).
- High Pressure Triaxial Automated System (HPETAS) which is based on High Pressure Controllers (≥ 16MPa).

Load Frames & Triaxial Cells

- 50kN and 100kN velocity controlled devices with serial PC connectivity.
- 100kN, 250kN, 400kN, 500kN, 1000kN velocity, position and direct load feedback control with serial or IEEE PC connectivity.
- 2000kPa, specimens up to 50, 100 or 150mm (load frames
 > 50kN for 150mm cell due to size).
- 3400kPa, specimens up to 77mm.
- 14MPa, specimens up to 38mm or 100mm.
- 20MPa, specimens up to 70mm.
- 64MPa, specimens up to 38mm or 100mm.
- 128MPa, specimens up to 50mm.



Fig 2. ETAS Temperature Controlled Features

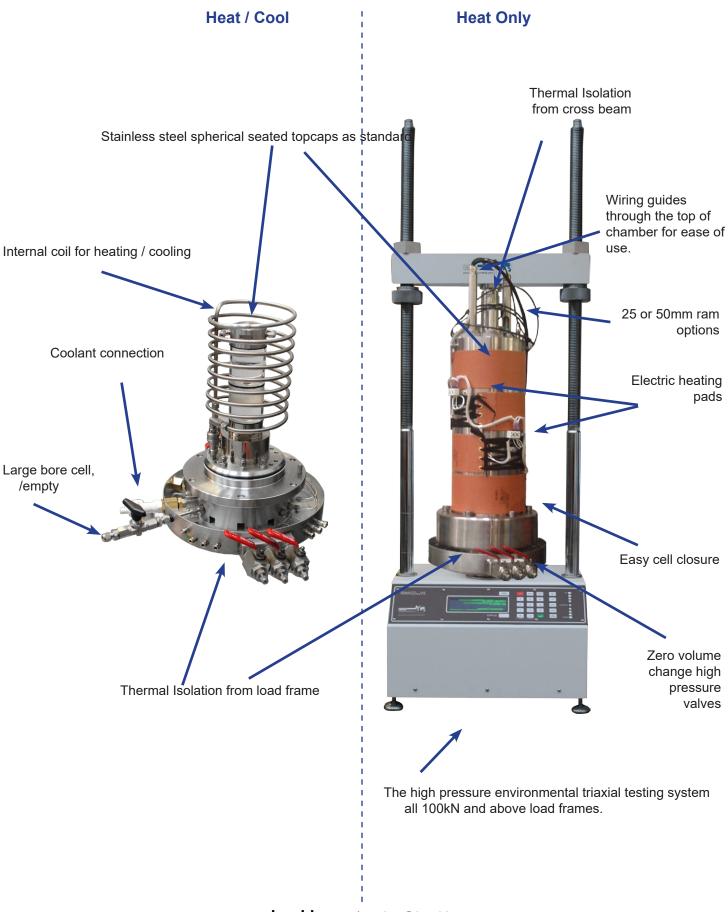

www.hoskin.ca | sales@hoskin.ca Vancouver | Oakville | Montréal



GDSLAB Control Software


GDSLAB is the control and data acquisition software for geotechnical laboratory applications. GDSLAB starts with a core application known as the kernel. The GDSLAB kernel allows for data acquisition from your hardware, but no test control. Simply add the appropriate module or modules to complete the test suite functionality you require. GDSLAB is compatible with all existing GDS equipment and furthermore key hardware from other manufacturers.

GDSLAB has the ability to be to your hardware of choice, no matter how unique the arrangement. A text (*.ini) or initialisation is created that describes the hardware connectivity to the PC. The hardware layout is available in graphical format via the GDSLAB 'object display'. This makes setting up the devices and checking the connectivity extremely simple.

GDSLAB REPORTS Presentation Software

GDSLAB REPORTS software presents data obtained by GDSLAB to the National Standard, BS 1377:1990. The program can be used to present data whether obtained from a GDSLAB data or inputted by hand.

GDSLAB Reports can as be used with other manufacturer's dataloggers as well as all versions of GDS data logger. The results can be exported as a CSV into Microsoft Excel, allowing the user to customise the layout of results.

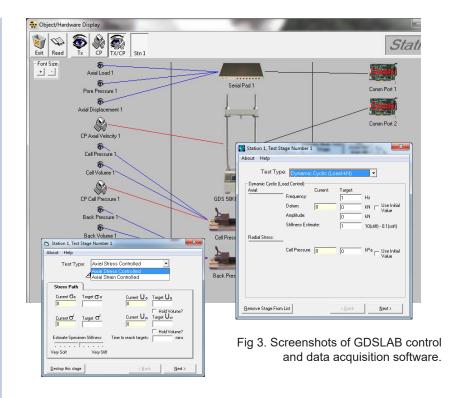

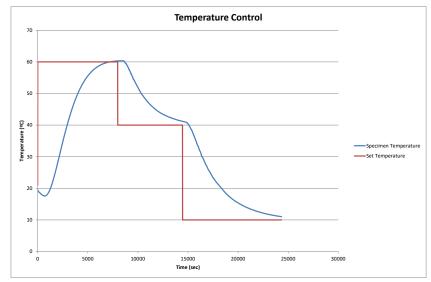


Fig 4. Example data for the rate of temperature change in a GDS HPETAS. Please note system will change the thermal properties and maximum rates of change of temperature.

Applications



Frozen Soils



Gas Hydrate

Images courtesy of : U.S. Geological Survey Department of the Interior/USGS U.S. Geological Survey/photo by Suzanne Weedman (Frozen Soils) & J. Pinkston and L. Stern (USGS - Gas Hydrate &)



Technical Specification

Load Frame & Triaxial Cell:

Load Range (kN): 50, 100, 250, 400, 500, 1000, 2000

Pressure Range (MPa): 8, 14, 16, 20, 32, 64, 100

Sample Sizes (mm): 38 to 150 diameter

Frame Computer Interface: USB

Temperature Control:

Combined Cooling / Heating: -20°C to +65°C, -30°C to +85°C, -30°C to +100°C

Heating Only System: Ambient to 65°C, Ambient to 100°C

Cooling Only System: Options to -30°C
Computer Interface: RS232 serial cable

Data Acquisition:

8 Channel Serial Pad:	May be used within any	of the system combinations.	This 16 bit device has 8 computer controlled
o Chamilei Seriai Pau.	gain ranges,	designed to suit transducer	s used in a triaxial test.

Pressure Controllers:

ELDPC:	Available in one	1MPa pressure and 200cc volume.
STDDPC:	Available in 200cc only, but w	ith pressure ranges from 100kPa to 4MPa (default is 3MPa).
ADVDPC:	Available in pressure ranges from 100kPa to 8000KPa and a volumetric capacity of 200cm³, or up to 2000kPa with avolumetric capacity of 1000cm³.	
HPDPC:		available in pressure ranges from of 8MPa to 150MPa and a ³ . High pressure controllers can be used with water, oil or air.

Software:

GDSLAB:	Control and acquisition software.
GDSLAB Reports:	A program to produce results from data obtained by GDSLAB to National Standards, BS 1377:1990, ASTM D2850, D4767, D7181.
Available Test Modules:	SatCon (saturation & consolidation), Standard Triaxial (strain-controlled shearing), K0 Consolidation, Permeability, Stress Paths (q-p and s-t linear paths), Advanced Loading, 4D Unsat (stress paths for unsaturated soil)

Upgrade Options:

Bender Elements:	The system can be upgraded to perform P and S wave bender element testing.
Unsaturated Testing:	The system may be upgraded to perform Unsaturated Testing using Method A or D (Method B on request).
Local Strain Transducers:	Hall Effect transducers may be used in water up to 1700kPa. LVDT transducers come in 2 versions low pressure (up to 3500 kPa) version for use in water & high pressure (up to 200 MPa) version for use in nonconducting oil.
Gaseous Back Pressures:	Air. Carbon Dioxide. Nitrogen & Methane.



Why Buy GDS?

GDS have supplied equipment to over 84% of the world's top 50 Universities:

GDS have supplied equipment to over 84% of the world's top 50 Universities who specialise in Civil & Structural Engineering, according to the "QS World University Ranking 2019" report.

GDS also work with many commercial laboratories including BGC Canada, Fugro, GEO, Geolabs, Geoteko, Golder Associates, Inpijn Blokpoel, Klonn Crippen, MEG Consulting, Multiconsult, Statens Vegvesen, NGI, Ramboll, Russell Geotechnical Innovations Ltd, SA Geolabs, SGS, Wiertsema and Partners to name a few.



Would you recommend GDS equipment to your colleague, friend or associate?

100% of our customers answered "YES"

Results from our post-delivery survey asked customers for feedback on their delivery, installation (if applicable), supporting documentation, apparatus and overall satisfaction with GDS. The survey ran for two years.



Made in the UK:

All GDS products are designed, manufactured and assembled in the UK at our in Hook. All products are quality assured before they are dispatched.

GDS are an ISO9001:2015 accredited company. The scope of this applies to the approved quality administration systems relating to the "Manufacture of Laboratory and Field Testing Equipment".



Extended Warranties:

All GDS apparatus are covered by a 12 month manufacturers warranty. In addition to the standard warranty, GDS comprehensive extended warranties for 12, 24 and 36 months, for peace of mind against any repairs in the future. The extended warranties can be purchased at any time during the 12 months of ownership.



GDS Training & Installation:

All installations & training are carried out by engineers. A GDS engineer is assigned to each order throughout the sales process. They will quality assure the apparatus prior to shipping, if installation has been purchased, install the apparatus on the customers site & provide the training.



Technical Support:

GDS understand the need for ongoing after sales support, so much so that they have their own dedicated customer support centre. Alongside their support centre GDS use a variety of additional support methods including remote PC support, product helpsheets, video tutorials, email and telephone support.

